

# SpectraAnalyzer

FOOD



**ZEUTEC**

EN






# What it does for you

The SpectraAnalyzer® food is the ideal solution for routine analysis of major quality parameters during food production.

In modern food processing operations, reliable and accurate analysis solutions are necessary to provide customers with products of highest and - what is most important – consistent quality. In order to be most competitive in the world market, consistent high yields, top quality and low production costs are the objectives that need to be achieved. Designed as a modular system, the SpectraAnalyzer® food solution presents the quantitative analytical results of major quality parameters (e.g. moisture, protein, fat, sugars, starch, acidity, salt, ash) and many others within 45 seconds. There is no need to manually condition the sample and extra reagents do not have to be used. This analyzer solution provides highly accurate quality control parameters at no extra cost. As a stand alone system the SpectraAnalyzer® food can be operated very easily and intuitively – even close to the production line. The rugged construction and unique optical sample/reference setup ensures reliable operation in environments with fluctuating temperatures, vibration and dust.

The SpectraAnalyzer® food analyzer comes with many ready to use calibrations and a powerful software package to facilitate calibration fine tuning and extensive and automated logging as well as database storage of the analytical results on the analyzer, within the company intranet and/or the internet.

## FEATURES INCLUDE

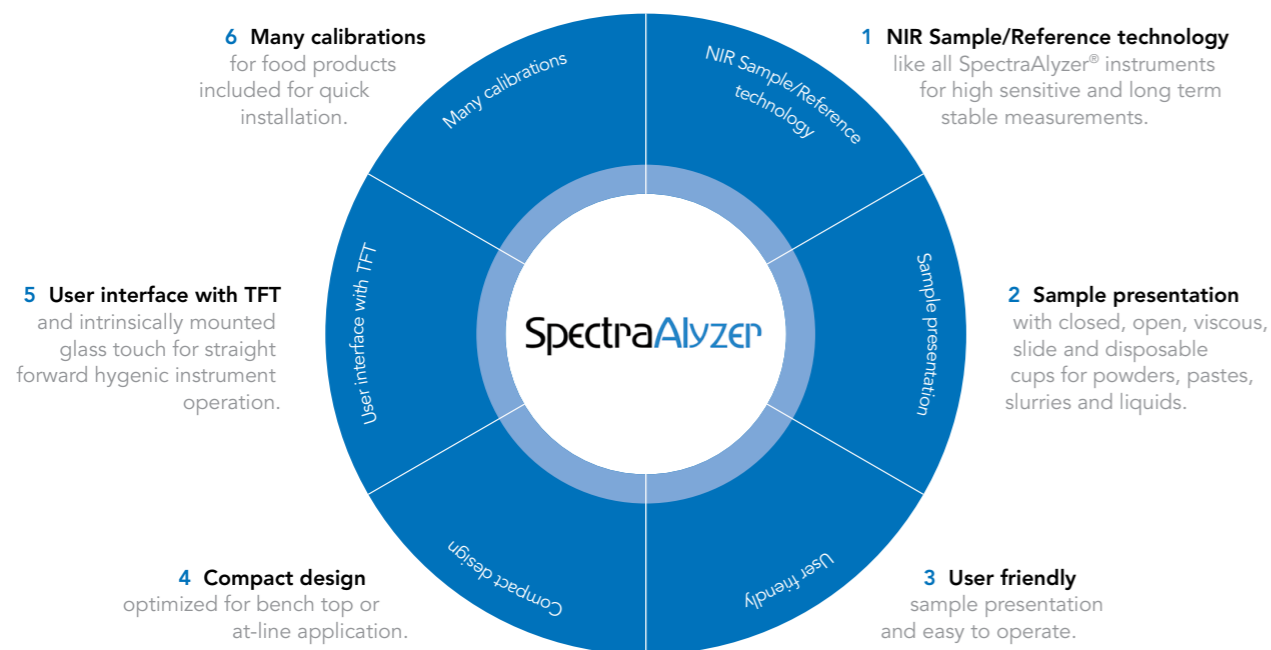

# Scope of supply

1 SpectraAnalyzer® food

2 Handbook, mains supply cable







# Technical data



## Design

For food products

Spectral range 1400 – 2400 nm

Dual beam system, Sample / Reference measurement

High signal to noise > 10000:1

Large expandable internal memory for calibrations, methods and history results

Auto-diagnostics

Graphical user interface, projected capacitive glass touch panel

## Accessories

Standard calibrations	Included
Keyboard	As optional accessory
Barcode reader	As optional accessory
Printer	Thermal ticket or office type printer as optional accessory
WiFi	Wireless LAN as optional accessory
AWX	Application worx, chemometric software package

## Specifications

Screen	TFT 800 x 480 pixel
Power requirements	min. 90 V AC (50 - 60 Hz) max. 260 V AC (50 - 60 Hz), 220 VA
Operating temperature	5 - 40 °C non-condensing
Interfaces	1 x front USB 2.0, 3 x USB 2.0, 2 x RS232, Ethernet, Audio out
Measurements	Height: 310 mm / Width: 300 mm / Depth: 480 mm
Weight	17 kg

## Reference

SpectraAlyzer® food	110-A100-1
---------------------	------------

### ZEUTEC Opto-Elektronik GmbH

Friedrich-Voß-Straße 11, 24768 Rendsburg

T (+49) 4331 - 136650

E [moreinfo@zeutec.de](mailto:moreinfo@zeutec.de)

[www.spectraalyzer.com](http://www.spectraalyzer.com)

# ZEUTEC

SpectraAlyzer® food Item No. 110-A100-1